

XPR Small Mass Comparators

Unrivalled Accuracy up to 10 kg



Unrivalled Accuracy

The innovative LevelMatic weighing pan eliminates corner load effects increasing the accuracy of your weighing results.




One-Click Certificate

The MC Link software automatically stores all your weight calibration data in a secure database. Calibration certificates can be issued with one click at the end of the weight calibration process.



Secure Data Transfer

The integrated USB and Ethernet interface makes it easy to connect to a PC or network. Data transfer to the mass calibration software MC Link is secure, fast and error-free.

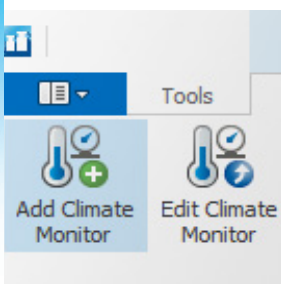

Effortless Compliant

Sophisticated monitoring of metrological relevant parameters within MC Link helps you ensure your results meet regulatory requirements.



Small Platform Comparators

Offer High Precision and Versatility

METTLER TOLEDO's XPR small platform mass comparators offer the highest weighing accuracy for determination of mass up to 10 kg. The combination of the innovative LevelMatic weighing pan and a dedicated draft shield delivers unmatched accuracy and repeatability. The state of the art user interface enables intuitive and efficient operation whilst the StatusLight™ built into the front of the terminal provides at-a-glance confirmation that it's safe to weigh and your results are valid. Our proprietary MC Link mass calibration software is supplied as standard. Enter your process specifications and you can be sure of secure results, full traceability and regulatory compliance.

XPR small mass comparators are the perfect choice for high-precision mass measurement applications, mass determination laboratories and mass calibration service providers.



XPR-SC Mass Comparators

	XPR2003SC	XPR2004SC	XPR5003SC	XPR10003SC
Material number	30317627	30317628	30317629	30317641

Limit values

Maximum capacity	2.3 kg	2.3 kg	5.1 kg	10.1 kg
Readability	1 mg	0.1 mg	1 mg	1 mg
Repeatability	1 mg	0.6 mg	2 mg	3.5 mg
Repeatability ABA (test load) ¹⁾	0.8 mg (2 kg)	0.25 mg (2 kg)	0.8 mg (5 kg)	1.5 mg (10 kg)
Linearity deviation	2 mg	1 mg	3 mg	7 mg

Typical values

Repeatability (test load)	0.7 mg (100 g)	0.3 mg (100 g)	1 mg (200 g)	2 mg (500 g)
Repeatability ABA (test load) ¹⁾	0.5 mg (2 kg)	0.15 mg (2 kg)	0.5 mg (5 kg)	0.8 mg (10 kg)
Linearity deviation	0.7 mg	0.5 mg	1.8 mg	3.2 mg
Eccentricity deviation (test load)	2.6 mg (1 kg)	Eliminated through LevelMatic		
Sensitivity offset (test load)	6.8 mg (2 kg)	6.8 mg (2 kg)	3 mg (5 kg)	34 mg (10 kg)
Settling time	3.5 s	3.5 s	3.5 s	3.5 s

Admissible Ambient Conditions

Temperature (°C)	10...30	10...30	10...30	10...30
Max. temperature change (°C/12h)	0.5	0.5	0.5	0.5
Relative humidity (%)	40 – 70	40 – 70	40 – 70	40 – 70

Dimensions

Weighing pan size (mm)	127 x 127	Ø 130	Ø 130	Ø 130
Object dimension (Ø, mm)	0...125	0...130	0...130	0...130
Object dimension (H, mm)	0...248	0...228	0...228	0...335

¹⁾ Repeatability calculated from 5 ABA cycles


Included Accessories

- MC Link Stand-alone mass calibration software (incl. 1 MT Instrument License)
- LevelMatic weighing pan (excl. XPR2003SC)

Selected Accessories

MC Link License – 1 MT Instrument	30208285
MC Link Upgrade – Multi-place License	30208289
MC Link Option – Barcode Package	30212767
MC Link Option – Audit Trail License	30208283
MC Link Validation Handbook	30212634
ClimaLog30 Certified	30078423
Reference Weight E1 2 kg Certified	159191
Reference Weight E1 5 kg Certified	159201
Reference Weight E1 10 Kg Certified	159211

Other Accessories on Request

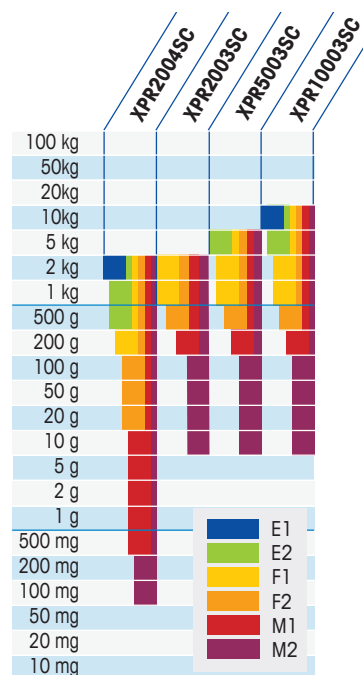

METTLER TOLEDO Group

Laboratory Weighing
Local contact: www.mt.com/contacts

www.mt.com/comparators

For more information

Subject to technical changes
© 06/2017 METTLER TOLEDO. All rights reserved
11795456A
Global MarCom 2302 RK

